

Anesthetic Gas Scavenging Inlet (Blower type)

SILBERMANN Anesthetic Gas Scavenging Inlet, blower type is designed to safely remove surplus anesthetic gas and explosive gas mixtures from the operating theatre.

The Inlet is a quick connect type, connected to a central scavenging system.



Anesthetic Gas Scavenging Blower type inlet and probe



Features

- Operation based on a central AGSS system.
- Internal adjusting valve ensures the required vacuum flow.
- Very quiet operation.
- Automatic shut off when probe is disconnected.
- Locking mechanism quick connect.
- Easy installation.
- Easy to maintain.

Applicable Standards:

- ISO 9170-2

Medical Gas System Standards:

- ISO 7396-2
- NFPA99
- HTM 02-01
- G0-1

Silbermann's Quality Management System:

- EN ISO 13485-2016
- ISO 9001-2015
- CE marked






Technical Data


- Outlet flow: 40 LPM
- Free air suction capacity: Adjustable 50 LPM
- Vacuum level: Adjustable 20-50 cmH₂O
- Exhaust pipe diameter: OD 7/8" nominal 20 cm long, K type.
- Material of construction:**
 - Box: Baked painted Aluminium
 - Front Panel: Aluminum



Dimensions for Recessed/Exposed [mm]

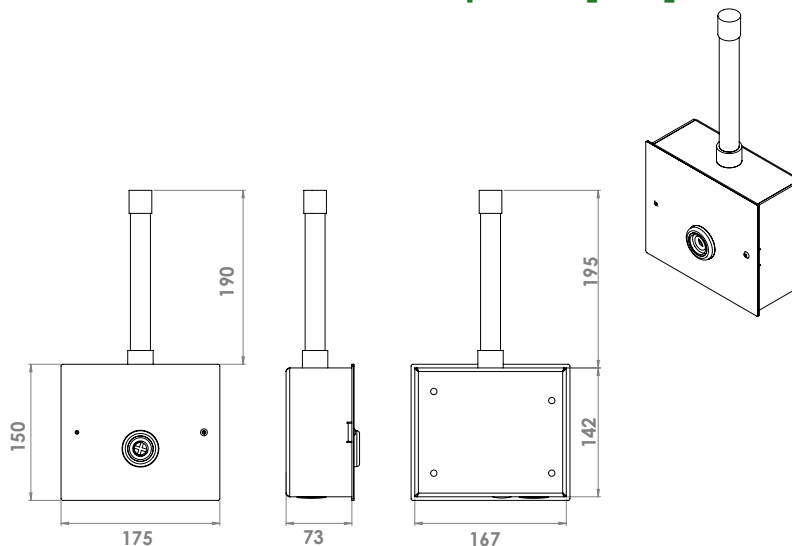

Dimensions for Pendant Mounted Cutout [mm]

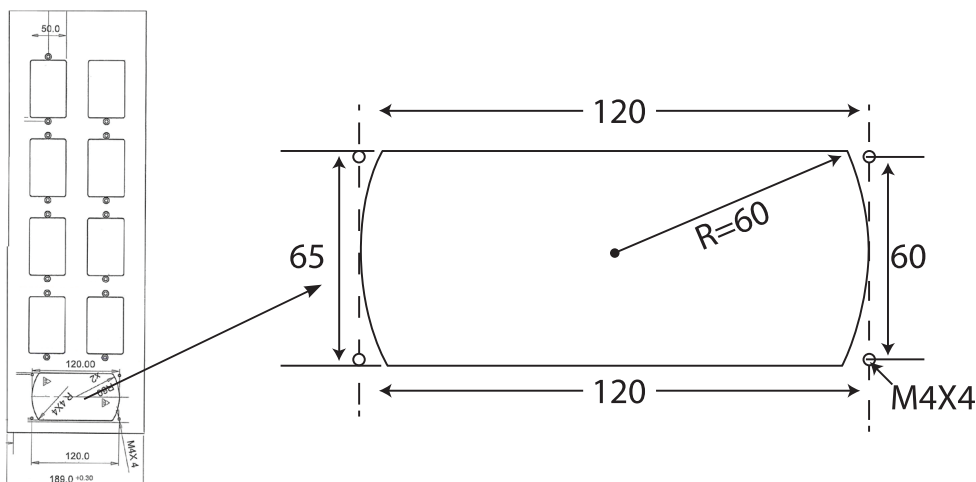

Anesthetic Gas Scavenging Inlet (Blower type)

Medical Gas Pipeline



Catalogue Numbers

Inlet Configuration	ISO 9170-2 Type
Recessed Box	3212880
Pendant/Panel Mounted	3212525



Related Products

Flexible Hose Assembly for OT

Hose Assembly 5 meter long includes:

- Probe
- Flexible hose - 5m long, 23 mm diameter.
- Rear adapter (to anaesthetic machine)

Type	Catalogue No.
Hose Assembly	2282200



Probe Only

Type	Catalogue No.
AGSS probe - EN	2282210

Distributed By:

SILBERMANN Technologies Ltd.

Address: #5 Harakevet st. POB 4605, Petah-Tikva, 4900715, ISRAEL

Tel: +972-3-9309186 | Fax: +972-3-9315721

E-Mail: info@silbermann.com, www.silbermann.com

